User's Manual
Manual del usuario

FreeStyle

Precision **H**

Blood Glucose and Ketone Monitoring System
Sistema de monitoreo de glucosa y cetonas en sangre

Definition of Symbols

These symbols are associated with your meter.

| | | | |
|--|---|---|---|
|  | Batch code |  | The batteries in this product should be removed and disposed in accordance with local regulations for separate collection of spent batteries. |
|  | Catalog number | | |
|  | Caution | | |
|  | Consult instructions for use |  | This packaging is capable of being recycled |
|  | Date of manufacture |  | Manufacturer |
|  | Do not drink |  | Sterilized using irradiation (lancets only) |
|  | Do not reuse |  | Temperature limitation |
|  | <i>In vitro</i> diagnostic medical device |  | Use by |
|  | Serial number | | |

Table of Contents

| | |
|---|----|
| Intended Use | 1 |
| Monitoring System Kit Contents | 2 |
| Setting Up the Meter | 3 |
| Enter Setup Mode | 3 |
| Beeper | 3 |
| Time | 3 |
| Date | 4 |
| Time Format | 4 |
| Date Format | 4 |
| Measurement Units | 5 |
| Preparing the Meter for Glucose Testing | 5 |
| Important Testing Notes | 5 |
| Calibration | 5 |
| Calibrating the Meter for Glucose Testing | 6 |
| Testing Blood Glucose | 9 |
| Important Testing Notes | 9 |
| Test Site | 9 |
| Performing a Blood Glucose Test | 10 |
| Understanding Blood Glucose Test Results | 13 |
| Low and High Blood Glucose Results | 13 |
| Preparing the Meter for Ketone Testing | 15 |
| Calibration | 15 |
| Calibrating the Meter for Ketone Testing | 15 |
| Testing Blood Ketone | 19 |
| Important Testing Notes | 19 |

| | |
|---|-----------|
| Test Site | 20 |
| Performing a Blood Ketone Test | 20 |
| Understanding Blood Ketone Test Results | 23 |
| High Blood Ketone Results | 23 |
| MediSense Glucose and Ketone Control Solutions | 25 |
| Important Control Solution Testing Notes | 26 |
| Performing a Control Solution Test | 27 |
| Managing Results | 30 |
| Viewing Test Results in Memory | 30 |
| Viewing Blood Glucose Averages | 33 |
| Transferring Results to a Computer | 35 |
| Error Messages | 35 |
| E-1 to E-9, What It Means, What To Do | 35 |
| Troubleshooting | 37 |
| Maintaining the Meter | 39 |
| Replacing the Battery | 39 |
| Cleaning the Meter | 41 |
| Specifications | 42 |
| References | 45 |
| Warranty Statement | 46 |

Intended Use

The FreeStyle Precision H monitoring system is intended for *in-vitro* diagnostic use by healthcare professionals. The system is indicated for the quantitative measurement of glucose (D-glucose) in fresh neonatal, venous, arterial, and fingertip capillary whole blood. The system is also indicated for the quantitative measurement of ketone (Beta-Hydroxybutyrate) in fresh capillary whole blood from the fingertip or venous whole blood samples provided the samples are used within 30 minutes of collection. The system is not intended for use in diagnosis or screening of diabetes mellitus, but is to be used as an aid in monitoring the effectiveness of diabetes control programs.

IMPORTANT: Use only **FreeStyle Precision H** Blood Glucose Test Strips and **FreeStyle Precision H** Blood β -Ketone Test Strips. Other test strips may produce inaccurate results.

See test strip instructions for use for more information about sample types.

Potential infection risk: Healthcare professionals performing blood tests with this system on multiple patients must always wear gloves and should follow the infection control policies and procedures approved by their facility.

Read the instructions in this user's manual. Failure to follow instructions will cause incorrect results. Practice the testing procedures before using the meter.

Observe caution when using around children. Small parts may constitute a choking hazard.

How FreeStyle Precision H Blood Glucose and Ketone Monitoring System Works

When you insert a test strip into the meter, the Apply Sample message shows on the meter's display window. When a blood sample or control solution sample is applied to the test strip, the glucose or ketone reacts with the chemicals on the test strip. This reaction produces a small electrical current that is measured. The result shows on the meter's display window.

IMPORTANT: Avoid getting dust, dirt, blood, control solution, water or any other substance in the meter's test strip port.

Monitoring System Kit Contents

Monitoring System Kit contains:

- | | |
|---|---|
| <ul style="list-style-type: none">• Meter | <ul style="list-style-type: none">• User's Manual |
|---|---|

Items not included:

- | | |
|--|--|
| <ul style="list-style-type: none">• Test Strips• Control Solutions• Data Management System | <ul style="list-style-type: none">• Lancing device• Lancets |
|--|--|

Setting Up the Meter

Enter Setup Mode

Start with the meter off (no test strip inserted).

Commands

- Press **ⓘ** to move between setup displays
- Press **▶** to move forward
- Press **◀** to go back
- Press **☀** to turn display light On and Off

Press and hold **ⓘ** until **✓** and **♪** appear. This indicates the beeper is turned on.

✗ and **♪** indicate the beeper is turned off.

Beeper

Press **◀** or **▶** to change the beeper setting.
Press **ⓘ** to save the beeper setting and move to set time.

The beeper sounds when:

- The test countdown starts
- The test countdown ends

Time

12:00^{AM}

The hour flashes. Press **◀** or **▶** until the correct hour appears. Press **ⓘ** to save the hour and to move to set minutes.

The minutes flash. Press **◀** or **▶** until the correct minutes appear. Press **ⓘ** to save the minutes and to move to set date.

Date

20 10

Year

The year flashes. Press ◀ or ▶ until the correct year appears. Press ⓘ to save the year and to move to set month.

10- 1

Month – Day

The month flashes. Press ◀ or ▶ until the correct month appears. Press ⓘ to save the month and to move to set day.

The day flashes. Press ◀ or ▶ until the correct day appears. Press ⓘ to save the day and to move to set time format.

Time Format

22:00

24 Hour Format

Note: If using the 12-hour time format, AM and PM appear.

The time format flashes. Press ◀ or ▶ to change the time format. Press ⓘ to save the time format and to move to set date format.

Date Format

10- 15

Month – Day

The date format flashes. Press ◀ or ▶ to change the date format. Press ⓘ to save the date format and to move to the measurement units.

15. 10

Day.Month

IMPORTANT: Please check that the correct time and date is set before you use the meter for the first time. This will help you keep records of when you monitor and will help make an informed decisions about your patient's care. **You must set the time and date to review averages.**

Measurement Units

mg/dL


The unit of measurement set by the factory appears. You cannot change this setting. Contact Customer Service for more information.

mmol/L

Press **1** to return to set beeper, or press and hold **1** to turn off the meter.

Preparing the Meter for Glucose Testing

Important Testing Notes

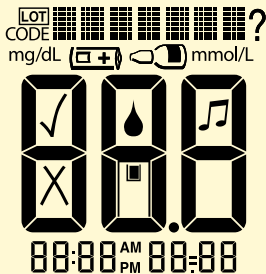
This display always appears when the meter is turned on. You should check that your meter matches the example you see here exactly every time your meter turns on.

Do not use the meter if the display check screen does not exactly match the example.

If the screen does not match, the meter may show an incorrect result. Please contact Customer Service.

IMPORTANT: This display means the meter is **not** calibrated for glucose testing.

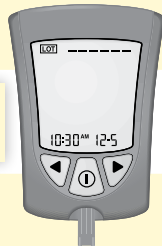Calibration

Calibrate the meter to:

- Match the test strips being used
- Obtain accurate results

Calibration is required:

- When using the meter for the first time
- When opening and using a new box of test strips


Calibrating the Meter for Glucose Testing

Remove the glucose calibrator (white) from the new box of test strips.

Peel the clear cover away from the corner with three bumps to open the calibrator package.

IMPORTANT: Use only the calibrator that is packaged in the box of test strips you are using.

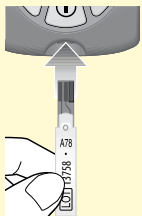Prepare the Meter

Hold the calibrator with the LOT number facing you. Insert the calibrator into the meter until it stops. This turns the meter on.

This display always appears when the meter is turned on.

Do not use the meter if the display check screen does not exactly match the example. Contact Customer Service.

See “Preparing the Meter for Glucose Testing” section for more information.

10:30^{AM} 12-5

Time and date (if set) appear next.

LOT 13758

The LOT number of the box of test strips you are using appears on the display. (Example shown)

The beeper (if sound is on) beeps.

Confirm the LOT Number

Check to see that the LOT number matches on all the following areas:

- Meter display
- Test strip calibrator
- Test strip instructions for use
- Test strip foil packet

IMPORTANT: Do not discard the calibrator until you have used all the test strips in the box.

| LOT Number | What It Means | What To Do |
|---|---|--|
| Matches all: LOT 13758 (Example) | The meter is calibrated for the box of test strips you are using. | You may now test your patient's blood glucose. |
| Does not match all | The meter may not be calibrated for the box of test strips you are using. | Check to see that you are using the calibrator that came packaged in the box of test strips you are using. Recalibrate the meter. |

IMPORTANT: Contact Customer Service if the LOT number still does not match. **Do not** test your patient's blood glucose. The meter may produce inaccurate results.

Remove the calibrator from the meter and store it in the box of test strips.

Recall the LOT Number

Press and continue to hold **▶** and **ⓘ** until all the following items appear on the display:

- Display check
- Glucose LOT number
- Time and date of the most recent glucose calibration

| LOT Number | What It Means | What To Do |
|-------------|---|--|
| Correct | The meter is calibrated for the box of test strips you are using. | Press and hold ⓘ for 2 seconds to turn off the meter. |
| Not Correct | The meter may not be calibrated for the box of test strips you are using. | Check to see that you are using the calibrator that came packaged in the box of test strips you are using. Recalibrate the meter. |

IMPORTANT: Improper calibration causes incorrect results.

Testing Blood Glucose

Important Testing Notes

E.g. Expiration date: March 31, 2012

- Use the meter within the test strip operating range
- Use the test strip immediately after removing it from its foil packet
- Use each test strip only once

Do not use expired test strips; they may produce inaccurate results. The expiration date is the last day of the month if only the month and year are shown.

Do not use a wet, bent, scratched or damaged test strip.

Do not use the test strip if its foil packet has a puncture or tear in it.

Note: See test strip instructions for use for operating range and more information on limitations of use.


Test Site


Fingertips

Note: Use fingertip blood samples for blood glucose testing; see test strip instructions for use for additional blood sample types.

Performing a Blood Glucose Test

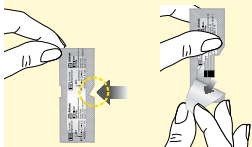Open the foil test strip packet at the notch and tear down to remove the test strip.

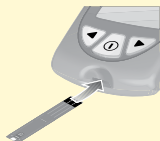Insert the Test Strip

With the contact bars (3 black lines) facing up, insert the test strip into the meter until it stops. This turns on the meter.

Note: The meter turns off after 3 minutes of inactivity. Remove and reinsert the unused test strip to restart the meter.

This display always appears when the meter is turned on.

Do not use the meter if the display check screen does not exactly match the example. Contact Customer Service.

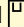
See “Preparing the Meter for Glucose Testing” section for more information.


10:30^{AM} 12-5

Time and date (if set) appear next.

LOT 13758


The LOT number of the box of test strips you are using appears next. The Apply Sample symbols  and  also appear, indicating the meter is ready for you to apply a sample to the test strip.

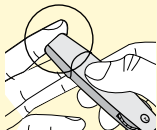Obtain a Blood Sample

Make sure that the sampling site is clean and dry before lancing. Use the lancing device to obtain a blood sample. (See lancing device instructions for use for more information.)

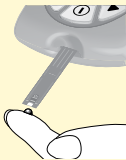Apply Blood to the Test Strip

Bring the blood drop to the white area at the end of the test strip. The blood is drawn into the test strip.

Hold the blood drop to the white area until the meter beeps (if sound is on) and the status bar appears on the display. This indicates the test strip has obtained enough blood.

:20

The countdown appears on the display as the meter checks the glucose level.

Note: Do not remove the test strip from the meter or disturb the test strip during the countdown.

If the countdown does not start, remove and discard the used test strip, turn off the meter and try again with a new strip.

mg/dL

104
9:49^{AM} 10-26

View the Result

The meter beeps (if sound is on) when the result appears on the display. The test is complete (example shown). The result is stored in memory.

Turn Off the Meter

Use the opened foil packet to remove the used test strip. This turns off the meter. Discard the used test strip properly.

Note: You can also press and hold **ⓘ** for at least 2 seconds to turn off the meter. The meter also turns off after 60 seconds of inactivity.

Understanding Blood Glucose Test Results

Low and High Blood Glucose Results

The meter displays results in mg/dL. The unit of measurement is preset. You cannot change this setting.

IMPORTANT: The meter displays results from 20–500 mg/dL. Low or high blood glucose results can indicate a potentially serious medical condition.

General:

Expected Results for Non-Diabetic, Non-Pregnant Adults

The expected glucose range for a non-diabetic, non-pregnant fasting adult is under 100 mg/dL. Two hours after meals, levels should be less than 140 mg/dL.⁷

Meter:

| Display | What It Means | What To Do |
|--|--|--|
|  Appears when result is: < 20 mg/dL | Severe low blood glucose or There may be a problem with the test strip | Repeat the test with a new test strip. |

Meter:

| Display | What It Means | What To Do |
|--|--|---|
| <p>KETONES?</p> <p>Appears with result which is: ≥ 300 mg/dL</p> | <p>Blood glucose level is high</p> | <p>Check blood ketone.</p> |
| <p>HI</p> <p>Appears when result is: > 500 mg/dL</p> | <p>Severe high blood glucose or There may be a problem with the test strip</p> | <p>Repeat the test with a new test strip. Check blood ketone.</p> |

Note: If you see the error message E-4, consult the Error Messages section in this user's manual.

IMPORTANT: If the blood glucose result appears to be inconsistent (lower or higher than expected), there may be a problem with the test strip. Repeat the test using a new test strip. Results that are incorrect may have serious medical consequences. Consult the prescribing physician before making any changes to diabetes medication plans.

Preparing the Meter for Ketone Testing

Calibration

Calibrate the meter to:

- Match the test strips being used
- Obtain accurate results

Calibration is required:

- When using the meter for the first time
- When opening and using a new box of test strips

IMPORTANT: This display means the meter is **not** calibrated for ketone testing.


Calibrating the Meter for Ketone Testing

Remove the ketone calibrator (purple) from the new box of test strips.

Peel the clear cover away from the corner with three bumps to open the calibrator package.

IMPORTANT: Use only the calibrator that is packaged in the box of test strips you are using.

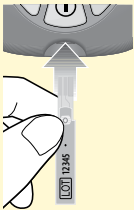Prepare the Meter

Hold the calibrator with the LOT/CODE number facing you. Insert the calibrator into the meter until it stops. This turns the meter on.

This display always appears when the meter is turned on.

Do not use the meter if the display check screen does not exactly match the example. Contact Customer Service.

See “Preparing the Meter for Glucose Testing” section for more information.

10:30^{AM} 12-5

Time and date (if set) appear next.

CODE 12345

The LOT/CODE number of the box of test strips you are using appears on the display. (Example shown)

The beeper (if sound is on) beeps.

Confirm the LOT/CODE Number

Check to see that the LOT/CODE number matches on all the following areas:

- Meter display
- Test strip calibrator
- Test strip instructions for use
- Test strip foil packet

IMPORTANT: Do not discard the calibrator until you have used all the test strips in the box.

| LOT/CODE Number | What It Means | What To Do |
|---|--|--|
| Matches all: CODE 12345 (Example) | The meter is calibrated for the box of test strips you are using | You may now test your patient's blood ketone. |
| Does not match all | The meter may not be calibrated for the box of test strips you are using | Check to see that you are using the calibrator that came packaged in the box of test strips you are using. Recalibrate the meter. |

IMPORTANT: Contact Customer Service if the LOT/CODE number still does not match. **Do not** test your patient's blood ketone. The meter may produce inaccurate results.

Remove the calibrator from the meter and store it in the box of test strips.

Recall the LOT/CODE number

Press and continue to hold ► and ⓘ until all the following items appear on the display:

- Display check
- Glucose LOT number
- Time and date of the most recent glucose calibration

After the last item appears, press ► or ◀ to display:

- Ketone LOT/CODE number
- Time and date of the most recent ketone calibration

| LOT/CODE Number | What It Means | What To Do |
|-----------------|--|--|
| Correct | The meter is calibrated for the box of test strips you are using | Press and hold ⓘ for 2 seconds to turn off the meter. |
| Not Correct | The meter may not be calibrated for the box of test strips you are using | Check to see that you are using the calibrator that came packaged in the box of test strips you are using. Recalibrate the meter. |

IMPORTANT: Improper calibration causes incorrect results.

Testing Blood Ketone

Important Testing Notes

E.g. Expiration date: March 31, 2012

- Use the meter within the test strip operating range
- Use the test strip immediately after removing it from the foil packet
- Use each test strip only once

Do not use expired test strips; they may produce inaccurate results. The expiration date is the last day of the month if only the month and year are shown.

Do not put urine on the test strip.

Do not use a wet, bent, scratched or damaged test strip.

Do not use the test strip if its foil packet has a puncture or tear in it.

Note: See test strip instructions for use for operating range and more information on limitations of use.


Test Site


Fingertips

Note: Use fingertip blood samples for blood ketone testing; see test strip instructions for use for additional blood sample types.

Performing a Blood Ketone Test

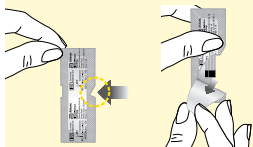Open the foil test strip packet at the notch and tear down to remove the test strip.

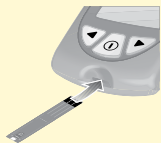Insert the Test Strip

With the contact bars (3 black lines) facing up, insert the test strip into the meter until it stops. This turns on the meter.

Note: The meter turns off after 3 minutes of inactivity. Remove and reinsert the unused test strip to restart the meter.

This display always appears when the meter is turned on.

Do not use the meter if the display check screen does not exactly match the example. Contact Customer Service.

See "Preparing the Meter for Glucose Testing" section for more information.

10:30^{AM} 12-5

Time and date (if set) appear next.


CODE 12345

The LOT/CODE number of the box of test strips you are using appears next.

KETONE

KETONE and the Apply Sample symbols  also appear, indicating the meter is ready for you to apply a sample to the test strip.

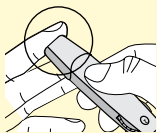Obtain a Blood Sample

Select a test site. Use the lancing device to obtain a blood sample. (See lancing device instructions for use for more information.)

Apply Blood to the Test Strip

Bring the blood drop to the white area at the end of the test strip. The blood is drawn into the test strip.

IMPORTANT: If the countdown does not start, you may not have applied enough blood to the test strip. See test strip instructions for use for re-application instructions. If the countdown still does not start, remove the used strip and discard it properly. Start a new test with a new test strip.

Hold the blood drop to the white area until the meter beeps (if sound is on) and the status bar appears on the display. This indicates the test strip has obtained enough blood.

KETONE

The countdown appears on the display as the meter checks the ketone level.

Note: Do not remove the test strip from the meter or disturb the test strip during the countdown.

: 10

View the Result

The meter beeps (if sound is on) when the result and **KETONE** appears on the display. The test is complete (example shown). The result is stored in memory.

Turn Off the Meter

Use the opened foil packet to remove the used test strip. This turns off the meter. Discard the used test strip properly.

Note: You can also press and hold ⓘ for at least 2 seconds to turn off the meter. The meter also turns off after 60 seconds of inactivity.

Understanding Blood Ketone Test Results

High Blood Ketone Results

The meter displays results in mmol/L. The unit of measurement is preset. You cannot change this setting.

IMPORTANT: The meter displays results from 0.0-8.0 mmol/L. Consult the prescribing physician before making any changes to diabetes medication plans.

Blood ketone is expected to be below 0.6 mmol/L.² High blood ketone may be caused by illness, fasting, vigorous exercise, or uncontrolled blood glucose levels.¹⁻³

Repeat a blood ketone test using a new blood ketone test strip when:

- HI appears on the display
- The result is unusually high
- You question the result
- The blood ketone result is 0.0 mmol/L, but the blood glucose is higher than 300 mg/dL

Meter:

| Display | What It Means | What To Do |
|--|---|--|
| Result is between: 0.6 – 1.5 mmol/L and the blood glucose result is: > 300 mg/dL | High blood ketone A problem requiring medical assistance may be occurring | Contact the prescribing physician for sick-day management instructions. |
| Result remains: 0.6 – 1.5 mmol/L or becomes: > 1.5 mmol/L | The patient may be at risk of developing diabetic ketoacidosis (DKA) ²⁻⁶ | Contact the prescribing physician immediately . |

Meter:

| Display | What It Means | What To Do |
|--|--|--|
| HI Appears when result is: > 8.0 mmol/L | Very high blood ketone or There may be a problem with the test strip | Repeat the test with a new test strip. If the result is HI, contact the prescribing physician immediately . |

MediSense Glucose and Ketone Control Solutions

IMPORTANT: When you open a control solution bottle for the first time, count forward 90 days and write this date on the control solution bottle using a pen that won't smear or wipe off. Throw away any remaining solution after this date.

Control solutions may be used with either glucose or ketone test strips.
Use control solution:

- To practice testing without using blood
 - To confirm that the meter and test strips work together properly
 - When you are unsure of your patient's blood glucose or blood ketone test results
- Contact Customer Service for information on how to obtain control solutions.

Important Control Solution Testing Notes

- Use only MediSense Glucose and Ketone control solutions with the meter
- Replace the cap securely on the bottle immediately after use

Do not use control solution past the expiration date.

Do not add water or other liquid to control solution.

Do not swallow or inject the control solution.


Do not use the control solution as eye drops.

IMPORTANT: Control solution results should fall within the control solution range printed on the test strip instructions for use. Check that the lot number printed on the test strip packet and instructions for use match.


- Repeat the test if control solution results are outside this range
- Stop using the meter if control solution results are continually outside the range printed on the test strip instructions for use. Contact Customer Service.

IMPORTANT: Control solution results do not reflect the blood glucose or blood ketone level.

ENGLISH Performing a Control Solution Test

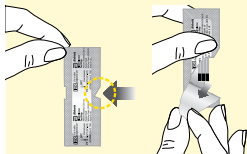Open the foil test strip packet at the notch and tear down to remove the test strip.

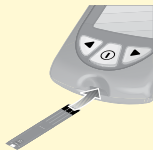Insert the Test Strip

With the contact bars (3 black lines) facing up, insert the test strip into the meter until it stops. This turns on the meter.

Note: The meter turns off after 3 minutes of inactivity. Remove and reinsert the unused test strip to restart the meter.

This display always appears when the meter is turned on.

Do not use the meter if the display check screen does not exactly match the example. Contact Customer Service.

See "Preparing the Meter for Glucose Testing" section for more information.

10:30^{AM} 12-5


Time and date (if set) appear next.

 13750

The LOT number of the box of test strips you are using appears next. The Apply Sample symbols  and  also appear.

Note: KETONE and  will appear if performing a ketone control solution test.

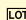 13750

Press  to mark the test as a control solution test.

The  appears. The meter is now ready for you to apply control solution to the test strip.

IMPORTANT: The test result will be saved to memory as a blood glucose or blood ketone result if not marked as a control solution test. This may affect the blood glucose averages.

Apply Control Solution to the Test Strip

Shake the control solution bottle to mix the solution. Apply a drop of control solution to the white area at the end of the test strip in the area shown. The control solution is drawn into the test strip.

Hold the control solution drop to the white area until the meter beeps (if sound is on) and the status bar appears on the display. This indicates the test strip has obtained enough control solution.

The countdown appears on the display as the meter checks the glucose or ketone level.

Note: Do not remove the test strip from the meter or disturb the test strip during the countdown.

:20

If the countdown does not start, remove and discard the used test strip, turn off the meter and try again with a new strip.

mg/dL

104

12:10 PM 10-20

View the Result

The meter beeps (if sound is on) when the result appears on the display. The test is complete (examples shown). The result is stored in memory as a control solution result.

Note: KETONE appears with the result if performing a ketone control solution test.

Compare the control solution result to the range printed on the blood glucose or blood ketone test strip instructions for use. The result should fall within the range.

Causes of out-of-range results include:

- The meter and test strips may not be working properly

Repeat the test with a new test strip. Follow the testing instructions carefully. Contact Customer Service if the result is still outside the printed range.

Turn Off the Meter

Use the opened foil packet to remove the used test strip. This turns off the meter. Discard the used test strip properly.

Note: You can also press and hold **1** for at least 2 seconds to turn off the meter. The meter also turns off after 60 seconds of inactivity.

Managing Results

Meter memory stores up to 450 events – control solution, blood glucose, and blood ketone results, and other meter information. The most recent events appear first.

Viewing Test Results in Memory

Start with the meter off (no test strip inserted).

Press **1** to turn the meter on.

This display always appears when the meter is turned on.

Do not use the meter if the display check screen does not exactly match the example. Contact Customer Service.

See "Preparing the Meter for Glucose Testing" section for more information.

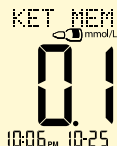
The most recent result, with time and date (if set) it was obtained, appears next. The following symbols indicate the type of result:

MEM = Blood glucose result

KET MEM = Blood ketone result

MEM  = Glucose control solution result

KET MEM  = Ketone control solution result

Press ◀ to review the previous result. Repeat this for each result you wish to review.

Note: To exit meter memory at any time, press and hold ⓘ.

MEM

Appear on the display when there are no previous results to view.

Press ▶ to return to a previously reviewed result. Repeat this for each result you wish to review.

Note: To exit meter memory at any time, press and hold ⓘ.

Viewing Blood Glucose Averages

IMPORTANT:

- Set the date and time before testing to be able to review averages
- Averages do not include glucose and ketone control solution results, blood ketone results, and results that do not show time and date
- Control solution results not marked as control solution tests may cause averages to be inaccurate
-  appear on the display when there are no current averages to view
-  blood glucose test results appear in averages as 20 mg/dL
-  blood glucose test results appear in averages as 500 mg/dL

While in memory, press . The 7-day average appears.

Note: n = number of results included in the 7-, 14-, or 30-day average.

Press ◀ to view the 14-day average.

Press ◀ to view the 30-day average.

Press ◀ to return to the 7-day average, or

Press ① to return to the memory screen.

Press and hold ① to turn off the meter. The meter turns off after 30 seconds of inactivity.

ENGLISH Transferring Results to a Computer

Transferring test results to a computer requires a compatible data management system. For more information, please go to our website or contact Customer Service.

Error Messages

| Message | What It Means | What To Do |
|---------|---|--|
| [- 1] | The temperature is too hot or too cold for the meter to work properly | <ol style="list-style-type: none">1. Move the meter and test strips to a location where the temperature is within the test strip operating range. (See test strip instructions for use for the appropriate range.)2. Wait for the meter and test strips to adjust to the new temperature.3. Repeat the test using a new test strip.4. If the error reappears, contact Customer Service. |
| [- 2] | Meter error | <ol style="list-style-type: none">1. Turn off the meter.2. Repeat the previous testing steps.3. If the error reappears, contact Customer Service. |
| [- 3] | There may be a problem with the test strip | <ol style="list-style-type: none">1. Review the testing instructions.2. Repeat the test using a new test strip.3. If the error reappears, contact Customer Service. |

| Message | What It Means | What To Do |
|---------|--|--|
| E-4 | The blood glucose level may be too high to be read by the system or There may be a problem with the test strip | <ol style="list-style-type: none">1. Repeat the test using a new test strip.2. If the error reappears, contact the prescribing physician immediately. |
| E-5 | Blood was applied to the test strip too soon | <ol style="list-style-type: none">1. Review the testing instructions.2. Repeat the test using a new test strip.3. If the error reappears, contact Customer Service. |
| E-6 | Calibration error or Test strip error | <ol style="list-style-type: none">1. Check the date setting on the meter.2. Check the expiration date on the test strip foil packet.3. Repeat the calibration using the calibrator that came with the test strip you are using.4. If the error reappears, contact Customer Service. |
| E-7 | Test strip may be damaged, used, or the meter does not recognize it | <ol style="list-style-type: none">1. Check that you are using the correct test strip for this meter. (See test strip instructions for use to verify the strip is compatible with this meter.)2. Repeat the test using a test strip for use with the meter.3. If the error reappears, contact Customer Service. |

| Message | What It Means | What To Do |
|------------------|---------------|---|
| E-8 or E-9 | Meter error | <ol style="list-style-type: none"> 1. Turn off the meter. 2. Repeat the previous testing steps. 3. If the error reappears, contact Customer Service. |

Troubleshooting

| | What It Means | What To Do |
|---|--|--|
| The meter does not enter test mode after inserting a test strip. | Test strip is not inserted correctly or fully into the meter | <ol style="list-style-type: none"> 1. With the contact bars (3 black lines) facing up, insert the test strip into the meter until it stops. This turns on the meter. 2. If the meter still does not enter test mode, contact Customer Service. |
| | No battery is installed Battery is installed incorrectly | <ol style="list-style-type: none"> 1. Install battery with (+) facing up. |
| | Dead battery | <ol style="list-style-type: none"> 1. Replace battery. Reset date and time, if necessary. |

| | What It Means | What To Do |
|---|--------------------------------------|--|
| The test does not start after applying the blood sample. | Blood sample is too small | <ol style="list-style-type: none">1. See test strip instructions for use for re-application instructions.2. If the countdown still does not start, remove the used strip and discard it properly.3. Repeat the test using a new test strip.4. If the test still does not start, contact Customer Service. |
| | Sample applied after meter turns off | <ol style="list-style-type: none">1. Remove the used strip and discard it properly.2. Review the testing instructions.3. Repeat the test using a new test strip.4. If the test still does not start, contact Customer Service. |
| | Defective meter or test strip | <ol style="list-style-type: none">1. Remove the used strip and discard it properly.2. Repeat the test using a new test strip.3. If the test still does not start, contact Customer Service. |


Maintaining the Meter

Replacing the Battery

The meter comes with a CR 2032 lithium (coin cell) battery installed. It provides power for about 1,000 tests.

- Appears on the display when the battery is low.
- Causes the display light to fail to turn on.

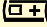When  appears, you may still use the meter and the results will be accurate. The display light will not work. Replace the battery immediately.

When  appears alone, the meter is not usable. The meter turns off automatically. Replace the battery immediately.

1. Gently push the battery cover in and up to remove it.

2. If the meter includes a plastic tab, pull on the tab to remove the battery.
If the meter does not include a plastic tab, push the battery upward to remove it.

3. Insert a new CR 2032, lithium (coin cell) battery: (+) facing up.

4. Slide the battery cover into place until it clicks.

Note: *Dispose of used batteries properly.*

5. Press and hold **I** to turn the meter on. If the meter does not turn on, check to see that the battery was installed properly.
If meter turns on, the meter may prompt you to reset the time and date. (See 'Setting Up the Meter' in this user's manual.)

Note: *Test results will not be lost even if the meter loses its time and date settings.*

Cleaning the Meter

Avoid getting dirt, dust, blood, control solution, or liquid in the meter test strip port.

Clean the outside of the meter using a damp cloth and mild soap.

Other acceptable cleaning solutions include: 10% bleach, 70% isopropyl alcohol or 10% ammonia

IMPORTANT:

Do not try to clean the test strip port.

Do not pour liquid into the test strip port or onto the buttons.

Do not immerse the meter in water or other liquid.

Specifications

| | |
|----------------------------|--|
| Size | |
| Length | 7.47 cm |
| Width | Top 5.33 cm |
| | Bottom 4.32 cm |
| | Depth 1.63 cm |
| Weight | 40 to 46 grams |
| Power Source | 1 CR 2032 Lithium (coin cell) battery |
| Battery Life | Approximately 1,000 tests |
| Memory | Up to 450 events, including control solution, blood glucose and blood ketone results, and other meter information. |
| Storage Temperature | |
| Meter | -13 ° to 131 °F (-25 ° to 55 °C) |
| Test strips | See test strip instructions for use. |

| | |
|--------------------------------------|--|
| Control solution | See control solution instructions for use. |
| Assay Method | Amperometric electrochemistry |
| Assay Range | |
| Blood glucose | See blood glucose test strip instructions for use. |
| Blood ketone | See blood ketone test strip instructions for use. |
| Control Solution Range | |
| Blood glucose | See blood glucose test strip instructions for use. |
| Blood ketone | See blood ketone test strip instructions for use. |
| Functions | |
| Blood glucose testing | |
| Blood ketone testing | |
| Memory | 450 events |
| Glucose averaging | 7-, 14-, and 30-day |
| Control solution marking and testing | |

| | |
|--------------------------|---|
| Beeper | On: When countdown starts and ends. |
| Display Light | On: Stays on during countdown; turns off 30 seconds after results appear. Stays on while transferring results to a computer. |
| Data Port | Yes |
| Operating Range | |
| Meter temperature | 50 ° to 122 °F (10 ° to 50 °C) |
| System temperature | See operating range of the test strip you are using. |
| | See test strip instructions for use. |
| Meter relative humidity | 10% to 90%, non-condensing |
| System relative humidity | 10% to 90%, non-condensing |
| System altitude | Clinical testing demonstrates that altitudes up to 7,200 feet (2,195 meters) above sea level do not affect results |

References

- ¹ Schade DS, Eaton RP. Metabolic and clinical significance of ketosis. *Special Topics in Endocrinology and Metabolism* 1982; 4:1–27.
- ² Wiggam MI, O’Kane MJ, Harper R, Atkinson AB, Hadden DR, Trimble ER, Bell PM. Treatment of diabetic ketoacidosis using normalization of blood 3-hydroxybutyrate concentration as the endpoint of emergency management. *Diabetes Care* 1997; 20:1347–1352.
- ³ Harano Y, Kosugi K, Hyosu T, Suzuki M, Hidaka H, Kashiwagi A, Uno S, Shigeta Y. Ketone bodies as markers for Type 1 (insulin-dependent) diabetes and their value in the monitoring of diabetes control. *Diabetologia* 1984; 26:343–348.
- ⁴ Ubukata E. Diurnal variation of blood ketone bodies in insulin-dependent diabetes mellitus and non-insulin-dependent diabetes mellitus patients: The relationship to serum C-peptide immunoreactivity and free insulin. *Ann Nutr Metab* 1990; 34:333–342.
- ⁵ Luzi L, Barrett EJ, Groop LC, Ferrannini E, DeFronzo RA. Metabolic effects of low-dose insulin therapy on glucose metabolism in diabetic ketoacidosis. *Diabetes* 1988; 37:1470–1477.
- ⁶ Hale PJ, Crase J, Natrass M. Metabolic effects of bicarbonate in the treatment of diabetic ketoacidosis. *Br Med J* 1984; 289: 1035–1038.
- ⁷ American Diabetes Association Position Statement: Diagnosis and classification of diabetes mellitus. *Diabetes Care* 2012;35(Suppl.1):S64–S71.

Limited Warranty

Abbott Diabetes Care Inc. (“Abbott”) warrants that the Abbott Diabetes Care Blood Glucose Monitor (“Meter”) shall be free from defects in material and workmanship for a period of two (2) years from the original date of purchase or delivery date, whichever is later, provided it is not modified, altered or misused. Under this Limited Warranty, if the Meter is defective in material or workmanship, Abbott’s sole obligation is to replace the Meter, free of charge, with a same or similar meter as determined by Abbott in its sole discretion.

Abbott warrants the performance of the Meter in accordance with its specifications if it is used as directed and provided that the failure to perform or misperformance of the Meter has not been caused in whole or in part by the use of test strips that are not the test strips manufactured by Abbott.

This Limited Warranty covers only the Meter and shall not apply to auxiliary equipment or disposable accessories, extends only to the original purchaser, and is not assignable or transferable.

Please see the User’s Manual for Instructions for Use, distributor and manufacturer information.

For warranty service, contact Customer Service for assistance and/or instructions for obtaining a replacement Meter. See User’s Manual for your Customer Service contact information. Abbott may require as a condition to obtaining warranty service that you return your Meter, postage prepaid, to an address specified by Customer Service.

TO THE EXTENT POSSIBLE UNDER LAW, THESE ARE ABBOTT’S COMPLETE EXPRESS WARRANTIES FOR THE METER AND STATE YOUR EXCLUSIVE REMEDIES, ABBOTT MAKES NO OTHER WARRANTIES, EXPRESS OR IMPLIED, AND ABBOTT EXCLUDES AND DISCLAIMS ANY OTHER WARRANTIES INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE AND ANY OTHER REMEDIES INCLUDING, BUT NOT LIMITED, TO COMPENSATIONS UNDER CONTRACT OR TORT LAWS. ABBOTT DOES NOT WARRANT THAT OPERATION OF THE METER WILL BE UNINTERRUPTED OR ERROR FREE.

Limitation of Liability

TO THE EXTENT POSSIBLE UNDER LAW, ABBOTT WILL NOT BE LIABLE FOR ANY INCIDENTAL OR CONSEQUENTIAL DAMAGES RESULTING DIRECTLY OR INDIRECTLY FROM USE OF THE METER OR FAILURE OF THE METER TO PERFORM IN ACCORDANCE WITH SPECIFICATIONS.

Some jurisdictions do not allow the exclusion or limitation of other express or implied warranties or incidental or consequential damages, so the above limitations or exclusions may not apply to you.

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Referencias

- ¹ Schade DS, Eaton RP. Metabolic and clinical significance of ketosis. *Special Topics in Endocrinology and Metabolism* 1982; 4:1–27.
- ² Wiggam MI, O’Kane MJ, Harper R, Atkinson AB, Hadden DR, Trimble ER, Bell PM. Treatment of diabetic ketoacidosis using normalization of blood 3-hydroxybutyrate concentration as the endpoint of emergency management. *Diabetes Care* 1997; 20:1347–1352.
- ³ Harano Y, Kosugi K, Hyosu T, Suzuki M, Hidaka H, Kashiwagi A, Uno S, Shigeta Y. Ketone bodies as markers for Type 1 (insulin-dependent) diabetes and their value in the monitoring of diabetes control. *Diabetologia* 1984; 26:343-348.
- ⁴ Ubukata E. Diurnal variation of blood ketone bodies in insulin-dependent diabetes mellitus and non-insulin-dependent diabetes mellitus patients: The relationship to serum C-peptide immunoreactivity and free insulin. *Ann Nutr Metab* 1990; 34:333–342.
- ⁵ Luzi L, Barrett EJ, Groop LC, Ferrannini E, DeFronzo RA. Metabolic effects of low-dose insulin therapy on glucose metabolism in diabetic ketoacidosis. *Diabetes* 1988; 37:1470–1477.
- ⁶ Hale PJ, Crase J, Natrass M. Metabolic effects of bicarbonate in the treatment of diabetic ketoacidosis. *Br Med J* 1984; 289; 1035–1038.
- ⁷ American Diabetes Association Position Statement: Diagnosis and classification of diabetes mellitus. *Diabetes Care* 2012;35(Suppl.1):S64-S71.

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in various jurisdictions. ©2010-2019 Abbott
FreeStyle and other brand marks are trademarks of Abbott Diabetes Care Inc.