

# Single-use Surgical Mask

SOFT & EASY TO BREATHE

Non-sterile




## Application

Worn by personnel to cover the mouth, nose and chin, the mask provides a physical barrier to block exhalation or ejection of pollutants from mouth and nose, and to reduce breathing in airborne particles.



## Product Introduction

Single-use Surgical Mask consists of three layers of nonwoven material:

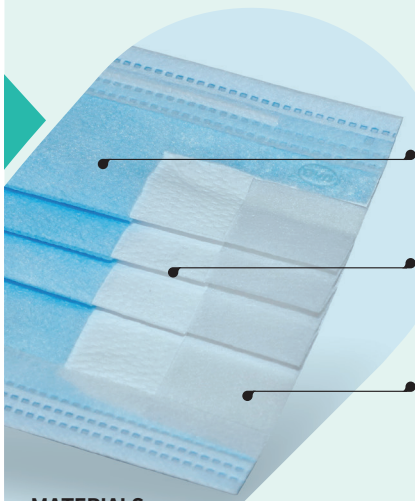

**Outer Layer**  
Polypropylene spunbond nonwoven fabric

**Middle Layer**  
Polypropylene melt-blown nonwoven

**Inner Layer**  
Polypropylene spunbond nonwoven fabric

### MATERIALS

Inner layer	Polypropylene spunbond nonwoven fabric
Middle layer	Polypropylene melt-blown nonwoven
Outer layer	Polypropylene spunbond nonwoven fabric
Ear loops	Polyester/nylon spandex blend
Nose clip	Metal wire with plastic covering



Non-sterile product. One-time use only.

## Wearing Instructions



1. Hold the mask by the ear loops with the colored side facing outwards and the nose bridge upwards, and place both ear loops around both ears.



2. Bend the nose clip to match the shape of the nose to prevent unfiltered air from entering.



3. Pull the mask to your chin to produce a tight seal.


# Single-use Surgical Mask

SOFT & EASY TO BREATHE

Non-sterile

## Product Specification

- 1 Mask length: 175 mm(6.89in.)
- 2 Mask expansion width: 165 mm(6.49in.). Post-stack width is 95 mm(3.74in.) with three stacks in the middle. Each stack width is not less than 10 mm(0.39in.). The nose clip is located on the upper edge of the mask, and the outside of the mask is darker in color.
- 3 The position of the upper nose clip is not more than 16 mm(0.63in.).
- 4 Nose clip: length > 90 mm(3.54in.); width approximately 3 mm(0.12in.)
- 5 Ear loop: 170 mm(6.69in.) in length and 3.5 mm(0.14in.) in diameter, composed of polyester, spandex and other materials, welded on the inner layer not more than 10 mm(0.39in.) from the edge.

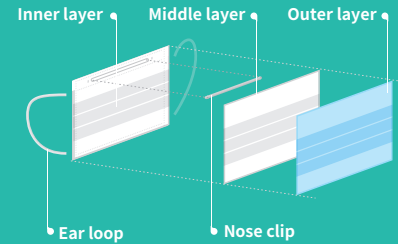

## Parameter of Single-use Surgical Mask

Name	Parameter
Product name	Single-use Surgical Mask
Material	Polypropylene spunbond nonwoven fabric, polypropylene melt-blown nonwoven, metal wire with plastic covering nose clip, polyester/nylon spandex blend ear loops
Product model	ASM2311X
Size	175 mm×95 mm / 6.89in.×3.74in.
Expiration date	2 years after production
Packaging specification	10pcs/bag, 50pcs/box, 2000pcs/case, 24 cartons/pallet
Storage	Store in a room with relative humidity level of no more than 80%, no corrosive gas and good ventilation to avoid high temperature.

## WARNING

- 1 The mask does not eliminate the risk of contracting any disease or infection.
- 2 Failure to properly use and maintain this product could result in illness or even death.
- 3 Only for one-time use, please destroy and dispose according to regulations after use.
- 4 Please use within the validity period.
- 5 DO NOT use if package is damaged.
- 6 Please refer to the instructions before use.
- 7 For adult use only.
- 8 Use this product immediately after the package is opened. Avoid hand contact within the interior part of the mask. Do not wash.
9. Intended for use in environments with moderate fluid exposure.
10. This face mask is authorized for distribution under FDA enforcement discretion only for the duration of the declaration that circumstances exist justifying the authorization of the emergency use of face masks under section 564(b)(1) of the Act, 21 U.S.C. § 360bbb-3(b)(1), unless the authorization is terminated or revoked sooner.

### Time use limitation

If the mask becomes damaged, soiled, or breathing becomes difficult, leave the contaminated area and replace and refit with a new mask.

### Important notice

To the extent permitted by law, Amsino shall not be liable for any loss or damage including any loss of business, loss of profits, or for any indirect, special, incidental or consequential loss or damage arising from reliance upon any information herein provided by Amsino. Nothing in this statement will be deemed to exclude or restrict Amsino's liability for death or personal injury arising from its negligence.

All information presented in this document is based on the latest data available at the time of printing. Amsino reserves the right to make changes at any time without prior notice.

## Approval and Standard

ASTM F2100 Level III

YY 0469-2011 standard

EN 14683:2019 Type IIR 