

EVAC-U-SPLINT®

Adult Vacuum Mattress

VERSATILE & ADAPTABLE Full Body Immobilization



Provides continuous spinal support & patient comfort in one complete package

The EVAC-U-SPLINT® vacuum mattress provides fast, effective and comfortable immobilization for the entire patient. Molding to the specific contours of the patient's body without applying pressure, the vacuum mattress provides uniform support and reduces pressure point tenderness.

The EVAC-U-SPLINT mattresses has six large padded handles secured to the perimeter of the mattress which provide improved security when moving and transferring a patient and also allows for rescuers to share handles when handing off a patient in tight quarters. The EVAC-U-SPLINT is filled with thousands of polystyrene beads which act as a thermal insulator to help maintain the patient's core body temperature decreasing the potential for hypothermia.



Optional backpack carry case available



QUICK AND EASY

Molds to each patient's unique body contours. Patient can be transported in a supine, knees flexed or a 30° elevated torso position.



INNOVATIVE VALVE

Includes quick connect and disconnect MaxiValve™. It is fast, effective and easy to use, even with gloves on. Audible click ensures a positive connection, even in low light.



SUPERIOR IMMOBILIZATION

Distributes body weight evenly eliminating specific pressure points. Ideal for injuries of the spine, back, hips, pelvis or trauma patients with multiple fractures.



VALUE AND PERFORMANCE

Integrated slats provide maximum support and the simple rescue-type strap system insures and security of the patient even if transported on their side for airway maintenance.



DESIGNED TO LAST

Made from high quality, extremely durable materials that are easy-to-clean and withstand years of use. Chosen for use in both summer and winter Olympic Games.



X-RAY COMPATIBLE


X-ray compatible.

HARTWELL MEDICAL

EVAC-U-SPLINT®

Adult Vacuum Mattress

Evacuation time is less than 25 seconds using the durable, very efficient all-aluminum manual pump with foot stirrup. Includes dual-wall braided hose.

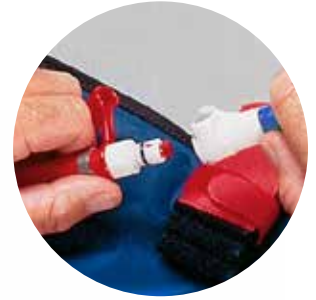

Features the simple, and intuitive MaxiValve™ that is extremely reliable and is guaranteed for the life of the product. Provides quick connect and disconnect. Includes protective leashed cap.



Features easy-to-apply color coded side-release buckles offering rapid and secure patient restraint.



SPECIFICATIONS

Weight:	
Mattress Only	12.1 lbs (5.5 kg) Model MT 1075-6
Mattress, Case, Pump	16.8 lbs. (7.6 kg) Model MT 3075-6
Dimensions:	
Unfolded	78.75"L x 29.5"W x 2.5"D (200 cm x 75 cm x 6.4 cm)
In Carry Case	31.5"L x 28"W x 6"D (80 cm x 71 cm x 15 cm)
Operating Temperature	-30°F to 150°F (-34°C to 66°C):
Load Capacity	350 lbs. (159 kg)
Fabric	Vinyl/polyester laminate
Filling	Polystyrene beads
Stiffener	Birch Plywood
Valve	Non-magnetic, self-sealing
Closures	Delrin side-release buckles
Carry Case	Polyester
Vacuum Pump	Aluminum cylinder with foot stirrup, stainless steel shaft

For optional ordering configurations please see our website.

*NSN-6515 01 408 7539



Contact your local authorized dealer or visit us online at www.HartwellMedical.com
For detailed product information, training and education materials please visit our website