

# GO<sub>2</sub>VENT

## GAS-OPERATED VENTILATOR



Hospital



• EMS/Transport




• MRI




No Electricity

### QUICK GUIDE

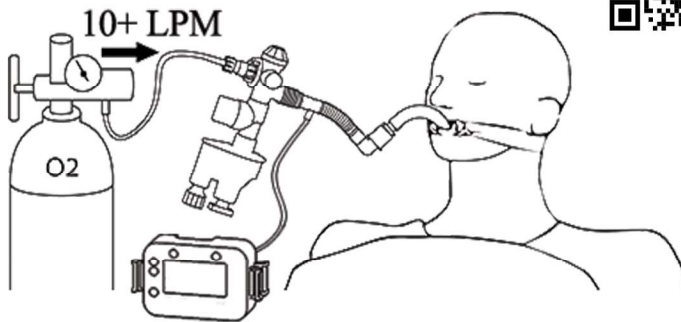

1. Connect GO<sub>2</sub>VENT and APM-Plus
2. Supply O<sub>2</sub> ≥10 LPM
3. Adjust FiO<sub>2</sub>, PIP / PEEP, Respiratory Rate (RR)

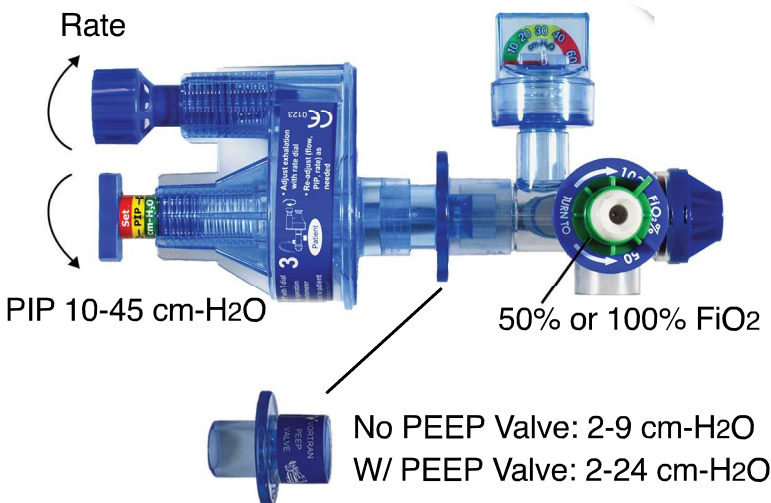

### A<sup>+</sup> APM-Plus

Airway Pressure Monitor

2 X AA



MONITOR:

- PIP
- PEEP
- RR
- I-Time
- E-Time
- I:E



ALARMS:

- High PIP
- Low PEEP
- High RR
- Non-Cycling

