

Security is always a priority...



Enhanced
version

Thomas Tube Holder


Whether you are securing an ET tube, LMA or other airway tube device, the Thomas Tube Holder is all you need! Designed to secure single and double lumen-airway tubes, the Thomas Tube Holder is fast and simple to apply and provides a reliable hold over tape fastening methods. In fact, the Thomas Tube Holder is ergonomically engineered to meet the needs of healthcare providers and offers a comfortable fit for both adult and pediatric patients.

The unique and enhanced design of the Thomas Tube Holder allows for the safe, effective security of ET Tubes, LMA's, Combitubes®, LTA's and other supraglottic airway tube devices. The Thomas Tube Holder gives the rescuer the confidence needed to secure the patient's airway... because security is always a priority.

Hook and loop padded strap provides security and comfort. Will not stretch or shrink when wet. Color-coded for ease of identification between adult and pediatric Thomas Tube Holder versions.

Larger access opening to allow for ease of suctioning and visual inspection of mouth.

Side stick facilitates fast, efficient application of the head strap around the patient's neck without movement of the head.

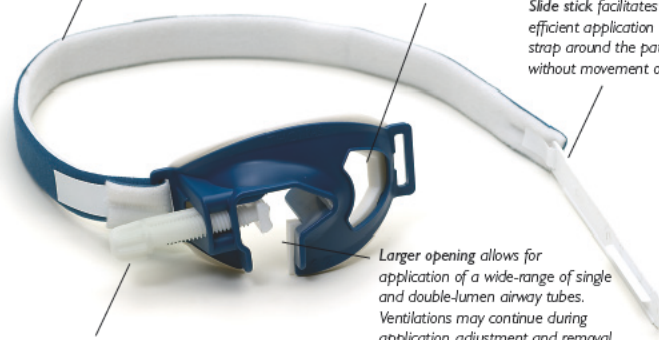


Larger opening allows for application of a wide-range of single and double-lumen airway tubes. Ventilations may continue during application, adjustment and removal of the Thomas Tube Holder.


Quick-set screw system securely holds tubes in place and facilitates fast, easy adjustment or removal of the Thomas Tube Holder.



Laerdal
helping save lives

www.laerdal.com

| Features | Benefits |
|--|--|
| Easy to Use | - Easy placement; eliminates the need for taping or tying the airway tube |
| Large opening  | - Accommodates single and double-lumen airway tubes including ET tubes, LMA's, Combitubes® and LTA's - Ventilations can continue during application, adjustment and removal of the Thomas Tube Holder |
| Large mouthpiece surface area | - Accommodates patients with removed dentures or facial trauma |
| Suction Port | - Large opening allows for view of patient's mouth and lip color - Easy access for suctioning and other oral care needs |
| Wide Head Strap | - Provides secure and comfortable fit; even when wet - Pink strap on pediatric Thomas Tube Holder for quick, easy identification in an emergency. |
| Pliable Bite Block | - Unique construction to prevent the patient from biting or pinching the tube and restricting airflow. Pliable material to protect the teeth and gums. |
| Versatile | - Pediatric and Adult sizes to accommodate the fastening of airway tube devices with outer diameters of between 4,3mm to 21mm |
| Convenient and Economic | - The Thomas Tube Holder is disposable – one time use - Prevents the need for messy clean-up and decontamination - Latex-free |

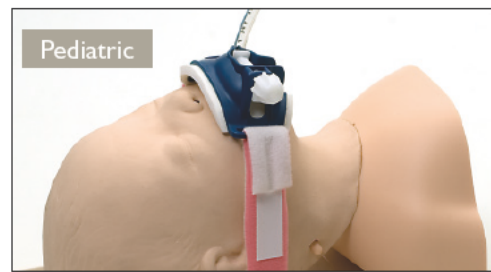

600-30000 Thomas Tube Holder Adult, Green (qty 1)

600-10000 Thomas Tube Holder Adult (qty 1)

600-20000 Thomas Tube Holder Pediatric (qty 1)



Laerdal

helping save lives