

Disposable intraosseous infusion needle



The disposable intraosseous infusion needle from CareFusion provides the ability to gain vascular access for delivering fluids and drugs to critically ill or injured children.

Intraosseous infusion (IO) is recommended only in life-threatening emergencies where vascular access is critical and conventional forms of access are not available or are unobtainable.



The disposable intraosseous infusion needle features:

- Sharp lancet point to facilitate entry into the bone
- Luer slip/luer lock connection for a secure fit to any syringe
- Adjustable depth guard to help control depth in needle insertion
- Two sizes to meet users' preference

Recommended needle placement:


- The proximal tibia may be used for children up to age 6
- The distal tibia is the preferred site for older children
- The distal femur

Ordering information

Catalog no.	Description	Qty/case
DIN1515X	15G (adjustable length $1\frac{5}{16}$ " to $1\frac{7}{8}$ " [24mm-48mm])	10
DIN1518X	18G (adjustable length $\frac{9}{16}$ " to $1\frac{7}{8}$ " [14mm-48mm])	10

CareFusion
Waukegan, IL

carefusion.com

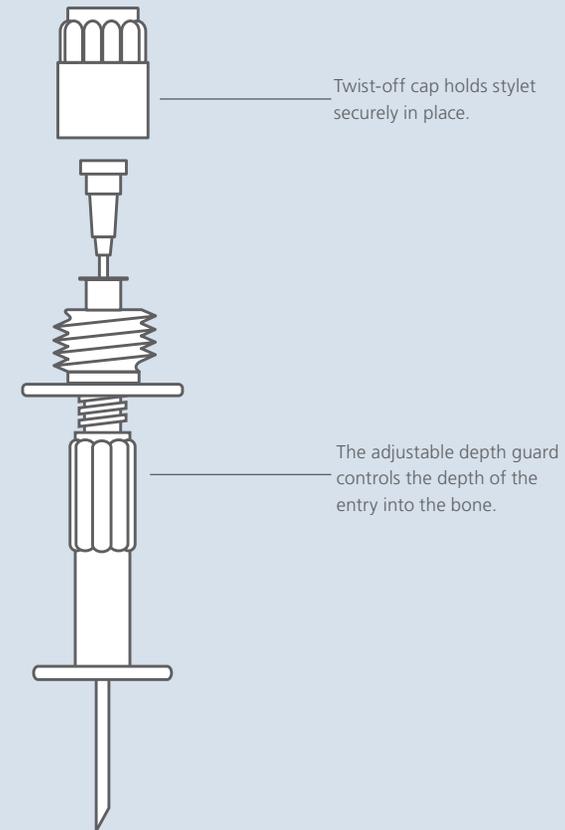

To place an order or for more information, please contact your Interventional Specialties Sales Representative or call 800.653.6827.

