



40 Parker
 Irvine, CA 92618
 (800) 326-4890
 Email: techservice-us@masimo.com

Masimo Accessory Warranty

Masimo warrants to the initial purchaser that each new SpO2 accessory and rainbow® cable will be free from defects in workmanship or materials for a period of six (6) months from the Lot Number date, plus one (1) month, to account for idle and shipping times. rainbow® sensors have a 1-year warranty.

- SpO2 Only: SpO2 only probes and RED LNCS & M-LNCS cables - 6 months
- rainbow®: Sensors, SpCO & SpMet – 1-year
- rainbow®: RC Cables – 6 months

Location of Alpha-Numeric Lot Number / Date Codes for Cables & Sensors:

- Label*: “DO NOT DISCARD” (attached to each cable & sensor):
- The lot number and Date Code may be stamped on the connector

MANUFACTURE YEAR: Use the first single-digit or two-digit number as the year (14, 15, 16, 17,18)

14= 2014 **15=** 2015 **16=** 2016 **17=** 2017 **18=** 2018


MANUFACTURE DATE CODES: Month (A-N)

A JANUARY **B** FEBRUARY **C** MARCH **D** APRIL **E** MAY **F** JUNE
G JULY **H** AUGUST **J** SEPTEMBER **K** OCTOBER **M** NOVEMBER **N** DECEMBER

Alpha date codes "I" and "L" are not used



Mfg. Date: 2007 - July

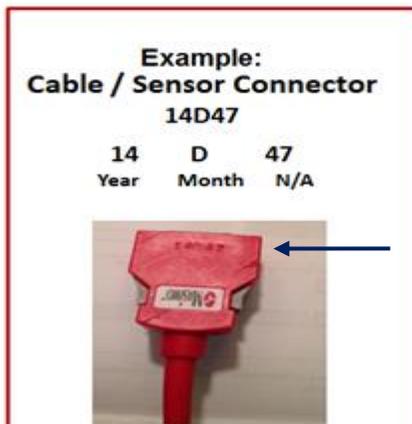

Mfg. Date: 2014 – April



Outer Box Lot Number

*Removal of the “Do Not Discard “
 Tag may void the warranty



40 Parker
Irvine, CA 92618
(800) 326-4890
Email: techservice-us@masimo.com

Contact Information:

- **Confirm the sensor or cable meets the 6 or 12 month warranty period**
- Contact Masimo Tech Support directly: (800) 326-4890 - Option (2) and Option (2) again.
- Masimo representative is available to perform sensor audits and recommendations at the local level.

rainbow® Part numbers covered by 1-year warranty:

- 2696- rainbow® Sensor DCI, SpO2/SpCO/SpMet 3 FT. 1/BX
- 2697- rainbow® Sensor DCIP, PED, SpO2/SpCO/SpMet 3 FT. 1/BX
- 3657- rainbow® Sensor DCI-6, ADT SpO2/SpCO/SpMet 6FT. 1/BX
- 3658- rainbow® Sensor DCI-P-6, PED SpO2/SpCO/SpMet 6FT. 1/BX

Direct-connect: **Switch to 2-Piece Solution: (RC-4 Cable #2406 + Sensor #2696 or #2697)**

- 2069- rainbow® Sensor DCIP-DC3, PED SpO2/SpCO/SpMet 3 FT. 1/BX
- 2202- rainbow® Sensor DCI-DC12, SpO2/SpCO/ SpMet 12 FT. 1/BX
- 2407- rainbow® Sensor DCI-dc8, SpO2/SpCO/SpMet 8 FT. 1/BX

rainbow® Part numbers covered by 6-month warranty:

- 2405- rainbow® RC-1 Ext. Cable 1-FT. 1/BOX
- 2406- rainbow® RC-4 Ext. Cable 4 FT. 1/BOX
- 3659- rainbow® RC-8 Ext. Cable 8 FT. 1/BX
- 2404- rainbow® RC-12 Ext. Cable 12 FT. 1/BOX, (Do not use)

Use the recommended solutions:

Harsh cleaning solutions will break down the cable jacket material.

CLEANING

1. Remove the sensor from the patient and disconnect it from the patient cable.
2. Clean the sensor by wiping with a 70% isopropyl alcohol pad or mild detergent.
3. Allow the sensor to dry thoroughly prior to placement on a patient.

or

1. If low-level disinfection is required, use a 1:10 bleach / water solution.
2. Saturate a cloth or gauze pad with the cleaning solution and wipe all surfaces of the sensor and cable.
3. Saturate another cloth or gauze pad with sterile or distilled water and wipe all surfaces of the sensor and cable.
4. Dry the sensor and cable with a clean cloth or dry gauze pad.

CAUTIONS:

- Do not use undiluted bleach (5% - 5.25% sodium hypochlorite) or any cleaning solution other than those recommended here because permanent damage to the sensor may occur.
 - Do not immerse the sensor or connector in any liquid solution.
 - Do not sterilize by irradiation, steam autoclave or ethylene oxide.
-