

NARC **iD**

NEW! Inventory Control System



CompX eLock®

NEW! NARC Box with NARC iD

The NARC iD Inventory Control System uses RFID tags and provides complete chain of custody information on narcotics, users, vehicles and consumption.

What does NARC iD do?

- ◆ Provides tracking information for narcotics
- ◆ Auto enrolling of drugs



Drug tote holds capsules in the correct position to be scanned by RFID tag on top of caps

Patent pending

PROVIDES DEA COMPLIANCE

CompX
SECURITY PRODUCTS

Life-Assist
Helping Heroes Save Lives

For more information, contact products@life-assist.com, call Life-Assist at 800-824-6016, or visit life-assist.com

Copyright 2019 © CompX Security Products

Any companies and/or products referred to herein are marks or registered trademarks of their respective companies, owners and/or mark holders.

LockView SOFTWARE

LockView 5iD Pro is easy to use, with many programmable options, and allows access to any eLock from any computer on the network.

LockView® software* features:

- ◆ Tracks narcotic consumed and ties it to the personnel accessing the NARC Box
- ◆ Automatically enrolls and translates the narcotic name into the software
- ◆ Manual inputs:
 - Drug manufacturer
 - Lot number
 - Expiration date
 - Date added
 - Date consumed



Minimum System Requirements

Windows Vista, 7, 8, 8.1, 10, Server 2008/2012 and R2 versions. 4GB RAM - Processor 2.0 GHz or higher. (8GB recommended) 4GB Hard drive space. SQL Server Standard, Enterprise, or SQL Express 2008/2012/2014 and R2 versions. If installing LockView on the Microsoft SQL Server please follow Microsoft SQL's minimum requirements. ODBC connection required to Server. Local Administrator Rights are required to install and register the software. .Net 4.0 Client Framework.

* LockView software sold separately

NARC iD features:


- ◆ Provides complete chain of custody information on narcotics, users, vehicles and consumption
- ◆ Provides DEA compliance
- ◆ Custom programmed tags provide automatic recognition of narcotic type and quantity – no user input required
- ◆ Tamper evident capsules hold narcotics and are color coded for easy identification
- ◆ Drug Tote ensures accurate reading of all RFID tracking tags
- ◆ Works with the 300 series CompX eLock
- ◆ Audit trails identify user, date, time of entry and narcotic consumed
- ◆ Allows input of the following narcotic conditions:
 - Incident number
 - Broken
 - Expired
 - Lost

Narcotics box features:

- ◆ Factory installed CompX eLock
- ◆ Heavy duty construction: seam welded 14 gauge steel
- ◆ Scratch-resistant powdercoated finish
- ◆ Built-in 12V-9V converter allows for dedicated power on-board EMS vehicles
- ◆ Auto-relocking door with mechanical key override
- ◆ Door position switch provides alert notifications

eLock features:

- ◆ Provides an audit trail of access attempts
- ◆ Uses existing credentials: Prox cards, Smart cards and Keypad access
- ◆ eLocks also available for retrofit on existing narcotics boxes and cabinets

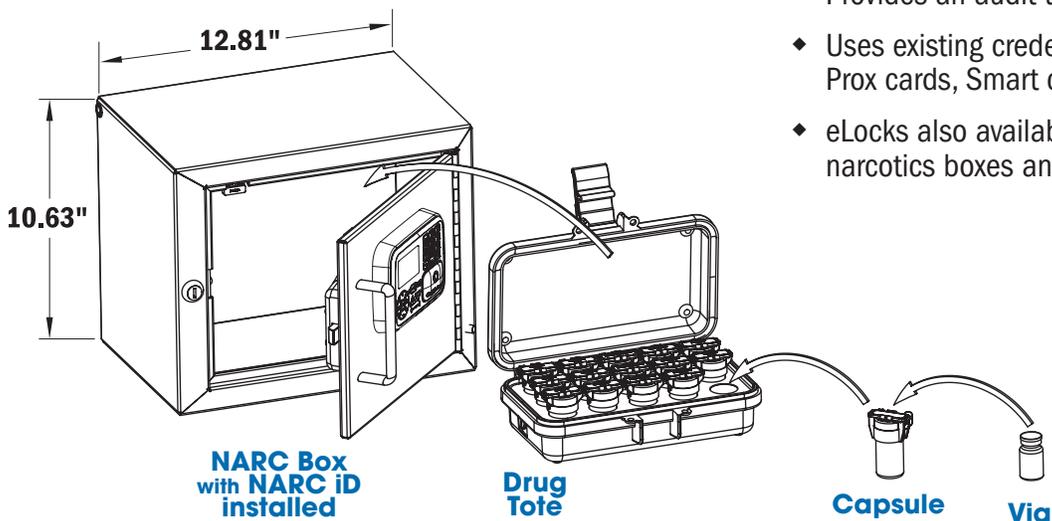

NARC iD features CompX eLock, an access control device that provides audit trail reporting and the ability to use existing card credentials.

Patent pending

CompX
SECURITY PRODUCTS



Reorder #: F2644-028

For more information, call **847.752.2500** or visit compxnarcid.com

Copyright 2019 © CompX Security Products / 847.752.2500 / compx.com / 715 Center St., Grayslake, IL 60030
Any companies and/or products referred to herein are marks or registered trademarks of their respective companies, owners and/or mark holders.