



Misty Max 10[™] Disposable nebulizer








Delivering fast and effective treatment, more efficiently. Accommodating all types of patients in all types of situations. These were the forces behind the development of the Misty Max 10™ nebulizer. Featuring a new type of jet design, the Misty Max 10™ nebulizer can speed up treatment and help ensure consistent medication delivery while offering you more flexibility in treatment options.

Nebulizer is available preassembled and fits virtually all common masks and adapters. All oxygen tubing is crush-resistant. Recommended for single patient use. The one-piece jet design and single-thread top connection make the unit easy to disassemble and reassemble.

Technical information (@ 8 L/min flow)*	
MMAD	1.61 µm
GSD	2.18
Respirable fraction % < 4.7µm (by mass)	85.2
Output rate	0.325 g/min
Treatment time	8 min, 9 sec
Residual volume	0.92 g
Fill capacity	10 mL

^{*} Please see the Misty Max 10 nebulizer technical data sheet for additional information, including performance characteristics at 5 L/min and 10 L/min flow.

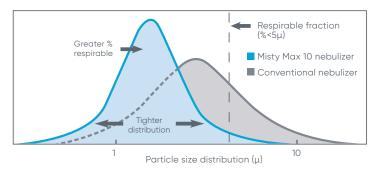

Particle size is an important characteristic of nebulizer performance. The Misty Max 10 nebulizer provides fine particles with a tighter distribution.*

* As proven by in-house testing

Ordering information

Cat. no.	Description
Nebulizers without mask, 50 per package	
002430	Nebulizer only
002431	With 7' oxygen tubing
002434	With baffled tee adapter, 7' oxygen tubing and mouthpiece
002435	With tee adapter, 7' oxygen tubing and mouthpiece
002438	With baffled tee adapter, 7' oxygen tubing, mouthpiece and 6" flextube
002446	With baffled tee, 7' U/Connect-It oxygen tubing, mouthpiece and 6" flextube
Nebulizer without mask, 25 per package	
002450	With patient "Y" mouthpiece, 7' oxygen tubing, 6" flextube and bacteria filter

Cat. no.	Description
Nebulizers with adult mask, 50 per package	
002433	With adult aerosol mask and 7' oxygen tubing
Nebulizers with pediatric mask, 50 per package	
002443	With pediatric aerosol mask
002444	With pediatric aerosol mask and 7' oxygen tubing
Masks only, 50 per case	
001261	Pediatric aerosol mask, short style
001263	Pediatric aerosol mask, under-the-chin style
001206	Adult aerosol mask

For more information, please contact your local Vyaire representative.

GLOBAL HEADQUARTERS

Vyaire Medical, Inc. 26125 N. Riverwoods Blvd. Mettawa, IL 60045 USA EC REP

EMERGO EUROPE Prinsessegracht 20 2514 AP The Hague The Netherlands

Vyaire Medical Pty Ltd Level 5, 7 Eden Park Drive Macquarie Park, NSW, 2113 Australia



Vyaire Medical, Inc. 26125 North Riverwoods Blvd Mettawa, IL 60045 USA



vyaire.com

For Global distribution.

© 2019 Vyaire Medical, Inc. or one of its affiliates. All rights reserved. Misty Max 10, Vyaire and the Vyaire Medical logo are trademarks or registered trademarks of Vyaire Medical, Inc or one of its affiliates. Medical devices class IIa according to Medical Devices Directive 93/42/EEC. Please read the complete Instructions For Use that come with the devices or follow the instructions on the product labelling. VYR-GBL-1900059

