

SSCOR DuCanto Catheter




“A simple and effective tool to support routine and emergency airway management.”

- James C. DuCanto, M.D.

- Larger internal diameter
- Designed to fit the Oropharyngeal Anatomy
- Can be used in the SALAD technique

The **SSCOR DuCanto Catheter™** is a single-use, non-sterile, rigid suction catheter specifically designed to be used by emergency care providers to aspirate liquids or semi-solids from a patient's upper airway.

It has a **larger inside diameter** than a standard Yankauer suction tip to facilitate the removal of fluids as well as solid material. The larger suction lumen makes it significantly less likely to clog in an emergency situation.

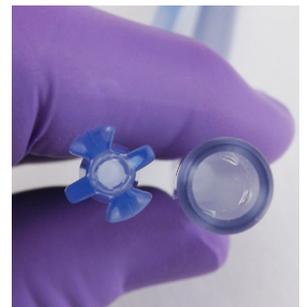

The **hyper curved catheter** is shaped more like a patient's airway making it easy to position the catheter and work with other complementary devices used to clear a patient's airway.



SSCOR DuCanto Catheter

Clearing The Airway Is Our #1 Priority

SSCOR.com

The **SSCOR DuCanto Catheter™** is constructed of a rigid PVC material. The larger internal diameter allows for larger pieces of emesis to pass while the lack of thumb port makes suction immediately available when the catheter is connected to the suction source. By removing the vacuum control, the caregiver does not need to locate and keep their thumb correctly positioned over the thumb port to clear a patient's airway.



SALAD

Suction Assisted Laryngoscopy
& Airway Decontamination

Instrumental in the **NEW SALAD Technique.**

The **SALAD (Suction Assisted Laryngoscopy and Airway Decontamination)** technique calls for constant suction provided by a catheter in the upper airway, continuously removing vomit, blood and other materials, around which the practitioner intubates the patient via either direct or video laryngoscopy.

Find more information on **SALAD**
at www.sscor.com/salad.

SSCOR DuCanto Reference/Part Number: 200-00002

Specifications:

Shaft diameter: .26" x .55" (6.6mm x 13.9mm)

Shaft length: 9.3" (236.22mm)

For best results, use **SSCOR #43200 OR #43203** patient connecting tubing. Both have a large internal diameter similar to the **SSCOR DuCanto Catheter**.

SSCOR | **DuCanto**
Catheter

Clearing The Airway Is Our #1 Priority



#43200 - 6 Feet



#43203 - 3 Feet

Dr. James C. DuCanto is an anesthesiologist with 20 years' experience in private practice who possesses substantial experience in the clinical practice and teaching of difficult airway management, as well as the use of patient care simulators.

Dr. DuCanto works at the Aurora St. Luke's Medical Center in Milwaukee, Wisconsin where he serves as a staff anesthesiologist, Director of the Simulation Center, and director of the Anesthesiology rotation for students and residents.

SSCOR, Inc.
11064 Randall St., Sun Valley, CA 91352
Email: sales@sscor.com or techsupport@sscor.com
Phone: (800) 434 5211 Domestic +1 818 504 4054 International
www.sscor.com & shop.sscor.com

CE 0413

SSCOR

Made in the USA of U.S. and Imported parts.
Lit 200-00002 Rev C