

Blizzard EMS Blanket (flat pack)

2-Layer Reflexcell® Protection Designed for EMS Use

- Compact, high-performing 2-layer Reflexcell® insulating blanket
- Regulates core body temperature to effectively prevent and treat hypothermia
- Protects from cold, wind, and rain
- Viable in temperatures as low as 18°F (-7.7°C)
- Compact, lightweight, waterproof, and vacuum-sealed packaging



Advanced Technology to Regulate Core Body Temperature

The Blizzard EMS Blanket is lightweight yet extremely effective. It's the only blanket of its kind designed to circulate heat across the entire body while retaining heat with unparalleled performance. This performance is achieved using Reflexcell® technology, a super-lightweight material that provides unprecedented thermal retention in the most demanding conditions.

Reflexcell is a framework of materials arranged to form cells that trap warm air and allow it to move throughout the entire blanket, ensuring that the patient is enveloped in warmth.

A network of elastic cords ensures that the blanket hugs the patient's body, stopping drafts and reducing convection, while a network of cells trap warm air, providing warmth and insulation where it is needed most—at the body's core. A metalized surface within these cells reflects heat back to the patient providing additional sustained warmth.

Blizzard EMS Blanket

(flat pack)


Product Highlights

Comprised of two layers of Reflexcell material, capable of withstanding blizzard-like conditions.

Prevents and treats hypothermia.

400% warmer than synthetic material blankets,
150% warmer than down-insulated materials.

Reflexcell® Technology

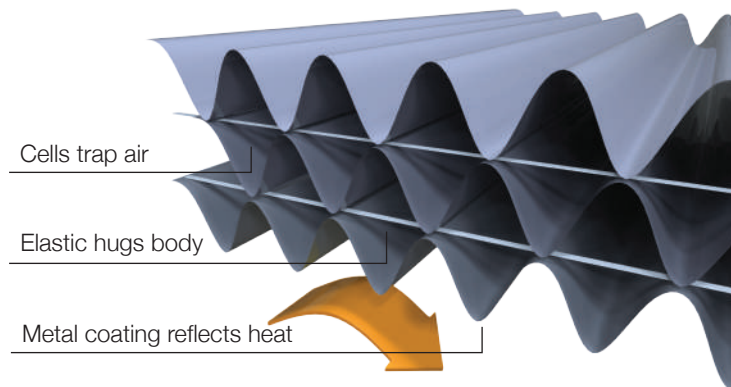

This compact, lightweight hypothermia blanket is an essential survival item of any emergency kit.

- Two layers of Reflexcell thermal retention material
- Elastic traps heat while metalized surface reflects heat back to user
- Compact, lightweight, waterproof packaging
- Sustains core body temperature to treat and prevent hypothermia
- Viable in temperatures as low as 18°F (-7.7°C)
- Seals with self-adhesive strip

Specifications

Description	2-layer Reflexcell insulating blanket, designed for EMS use
Part Number	BPS2-10FP
Packed Size	4.7in (12cm) wide 4.3in (11cm) long 1.9in (5cm) high
Size when open	7.6ft (2.33m) long 5.2ft (1.6m) wide, working 6.4ft (1.96m) wide, fully stretched
Weight	242g
Shelf Life	10 years
Colors	Orange

PerSys Medical
5310 Elm St.
Houston, TX 77081 USA
1.888.737.7978
info@ps-med.com

www.persysmedical.com



Life-Assist
Helping Heroes Save Lives