


NomoLine-0™ Capnography Sampling Lines

Airway Adapter Sets and Cannulas for Effective Capnography Monitoring

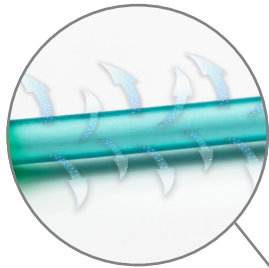


Patented moisture-wicking technology allows water evaporation from the sampling line

- > Water vapor-free sample provides accurate gas measurements, providing clinicians with confidence that EtCO₂ values reliably reflect the patient's physiological status



Single-twist luer connector for compatibility with many common monitors in numerous care areas

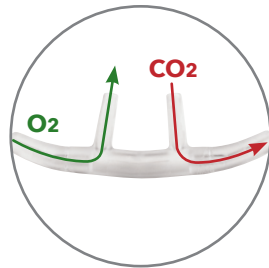

Split-flow cannula design reduces gas mixing, facilitating accurate EtCO₂ measurements

- > Cannula design delivers O₂ to the patient directly through the nasal prong



Ergonomically designed cannulas for patient comfort

- > Lightweight, soft cannulas with curved design fit the contours of the face to minimize pressure




NomoLine-O Features




Large, diverse product line allows use over all patient populations

- > For intubated patients, airway adapter sets with low dead space, allowing for fast rise time
- > For non-intubated patients, cannulas designed to facilitate EtCO₂ measurement in a variety of clinical care settings

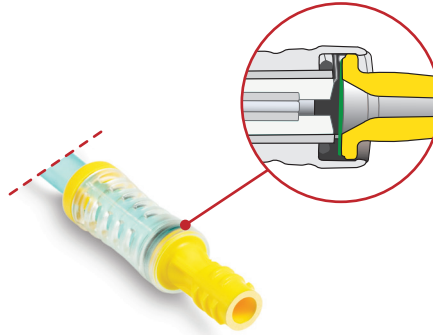

Hydrophobic bacterial filter reduces contamination of the gas analyzer

- > Extends the life of the module by preventing water intrusion

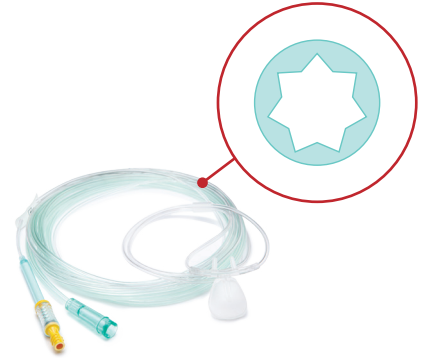

Star Lumen O₂ delivery tubing allows for continuous, accurate O₂ sampling and supply even if kinked

NomoLine-O Specifications

Length Options	.2 m and 4 m
Sampling Flow Rate	50 ml/min
Inline Bacterial Filter	BFE ≥ 99.9%
Tubing Inner Diameter	1 mm
Supplementary Oxygen Delivery	Up to 5 lpm
Airway Adapter dead space	
Adult/Pediatric	≤ 6 ml
Infant and Infant/Neonatal	≤ 0.7 ml
Material	Not made with natural rubber latex

Caution: Federal (USA) law restricts this device to sale by or on the order of a physician. See instructions for use for full prescribing information, including indications, contraindications, warnings, and precautions.

NomoLine-O is not licensed for sale in Canada.

Masimo U.S.
Tel: 1 877 4 Masimo
info-america@masimo.com

Masimo International
Tel: +41 32 720 1111
info-international@masimo.com

