


DISPOSABLE STRAPS

ECONOMY DISPOSABLE POLYPROPYLENE STRAPS

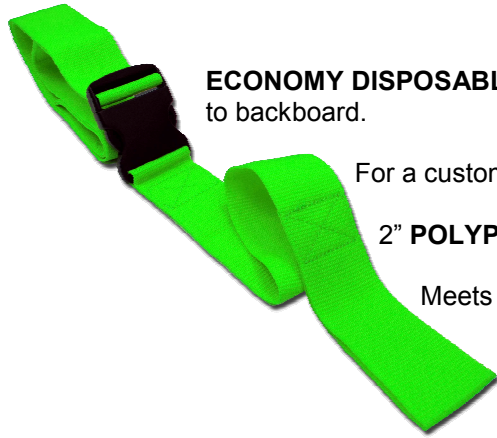


ECONOMY DISPOSABLE POLYPROPYLENE backboard straps are made to secure patients to backboard.

For a custom fit, economy disposable poly straps are easily adjustable.

2" **POLYPROPYLENE** webbing (1800 pound break strength)

Meets FDA requirements.

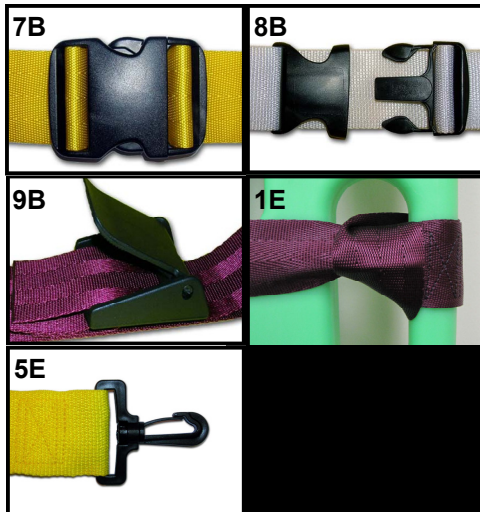

INSTRUCTIONS:

- 1.) Inspect strap before use.
- 2.) Position patient.
- 3.) Strap should have proper tension with no twists.
- 4.) Patient should be flat on backboard/stretcher.
- 5.) Wrap strap around patient and backboard/stretcher.
- 6.) Buckle and tighten strap to secure.

Intended as a guide only.

Consult your local medical agency for all protocols.

| | |
|--|--|
| 7B Plastic Side Release Buckle 400 Pound Break Strength | 1E Loop Lock Ends |
| 8B Plastic Double Adjust Buckle 400 Pound Break Strength | 5E Plastic Swivel Speed Clip Ends 140 Pound Break Strength |
| 9B Plastic Cam Buckle 5000 Pound Break Strength | |

| | | | |
|-------|--------|------------|------------|
| Black | Orange | Neon Green | Red |
| Blue | Green | Yellow | Olive Drab |
| White | | | |



For use by qualified and trained personnel only.